

Dosimetric evaluation of an efficient skin marker for marking scars, sarcomas, and larger treatment areas during proton therapy

Jenna Jatczak, MS, CMD; Thomas Houser, RT(T), CMD; Eric Kusmaul, RT(T) CMD; Shifeng Chen, PhD; Matthew Ferris, MD, Byongyong Yi, PhD; Nrusingh C. Biswal, PhD

Department of Radiation Oncology, University of Maryland School of Medicine, Baltimore, MD
Maryland Proton Treatment Center, Baltimore, MD

PURPOSE:

- A new skin marker, that helps to exclude skin from the automatic body contour is evaluated for the clinical use for proton beam therapy (PBT), for marking field borders, scars, sarcomas, and larger treatment areas.
- The dosimetric advantages and clinical efficiency are presented.

METHODS:

- An anthropomorphic phantom with the current marker on one side and the new marker on the other side was scanned with 3.0 mm slice thickness.
- The image set was duplicated, with marker densities overridden to air in one set and without override in the other set.
- Structures were contoured as per our clinical guidelines.
- Five monoenergetic plans, comprising the energy ranges in our clinic, with and without range shifters (70 -245 MeV) were generated on each image set, with a single beam at gantry angle 0° and snout position at 42 cm for no range shifter and snout at 16 cm for the 3 cm range shifter plans.
- For each plan, the field size was 20cmx20cm, 0.25 cm voxel size, 0.25 cm spot spacing, and 0.4 cm spot size with 10 MUs per spot.

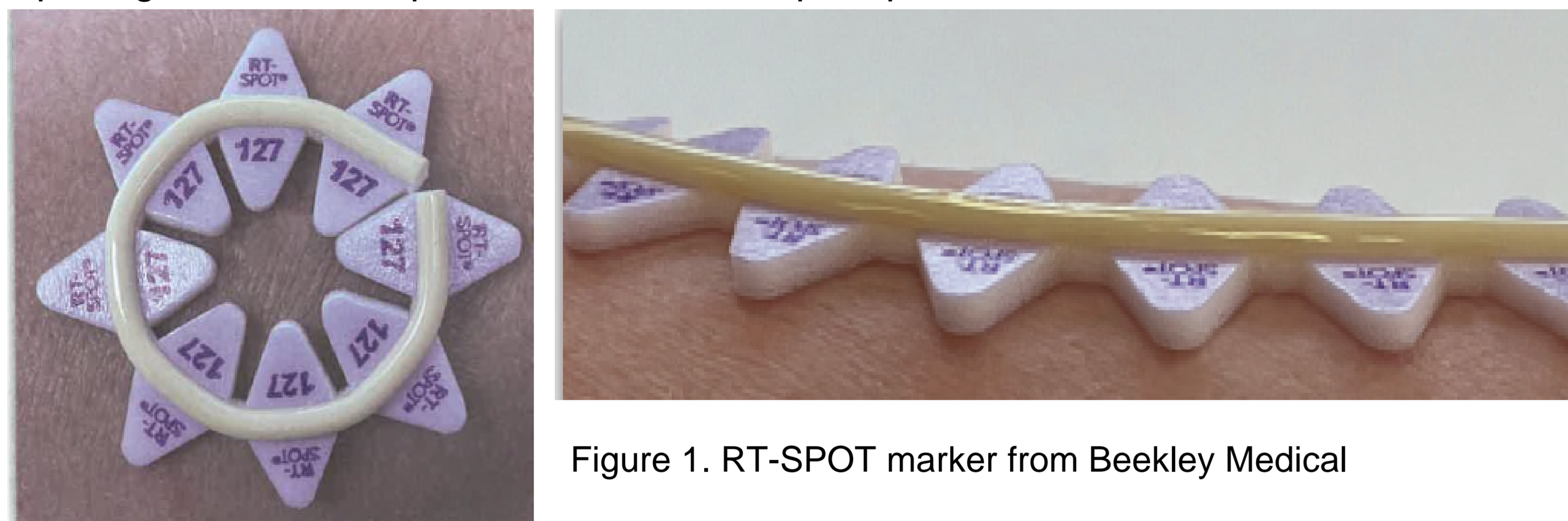


Figure 1. RT-SPOT marker from Beekley Medical

RESULTS:

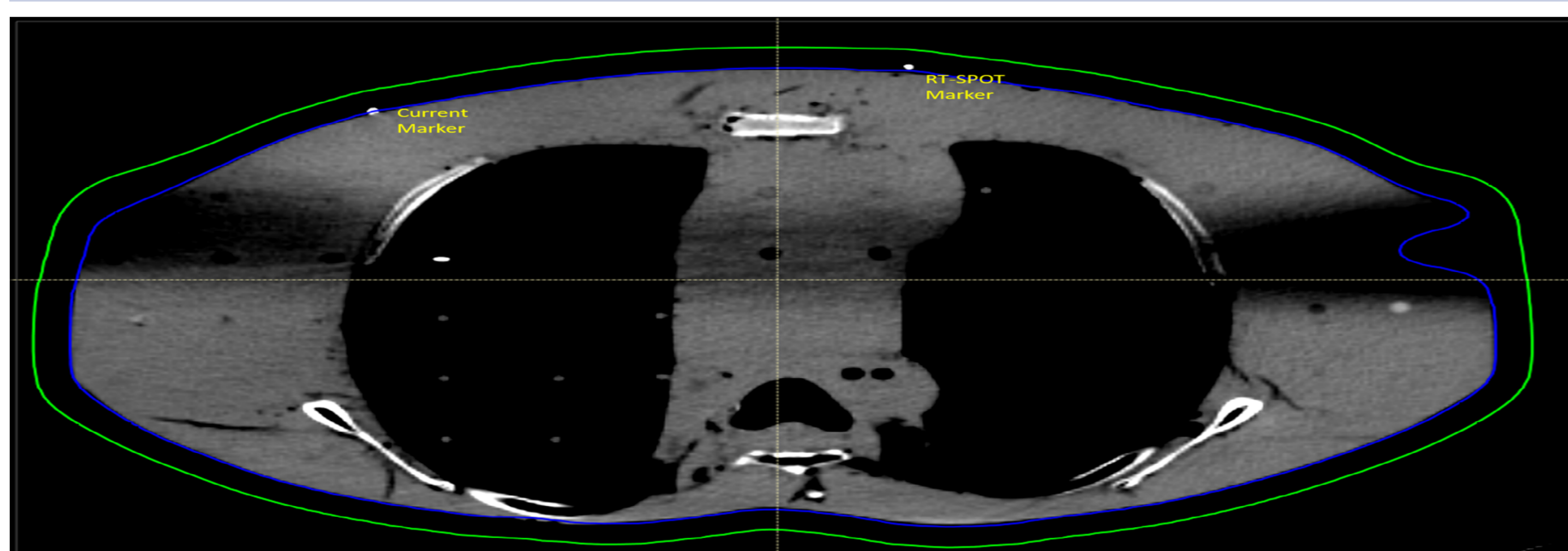


Figure 2. A representative slice of the CT image set, current marker on the right side and the RT-SPOT marker on the left side.

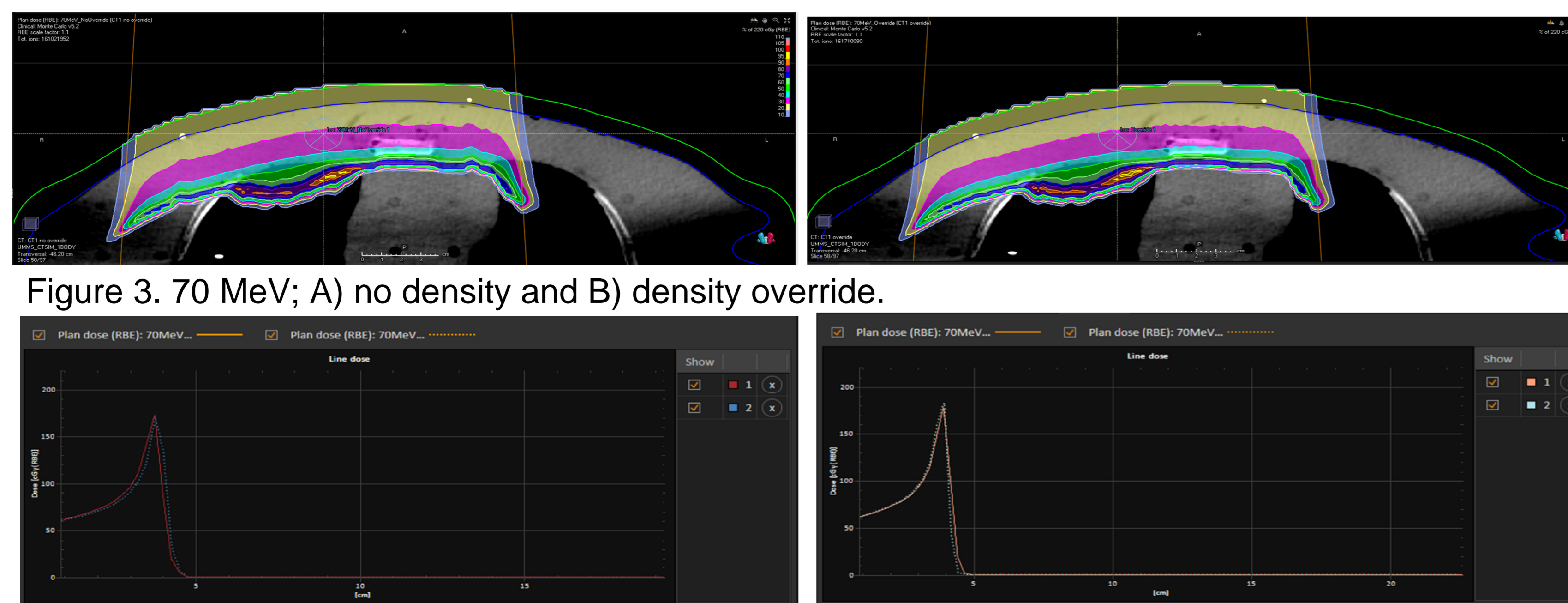


Figure 3. 70 MeV; A) no density and B) density override.

Figure 4. Dose profile perturbation at 70 MeV. A) current marker without (1) and with density override (2), B) RT-SPOT marker without (1) and with density override (2).

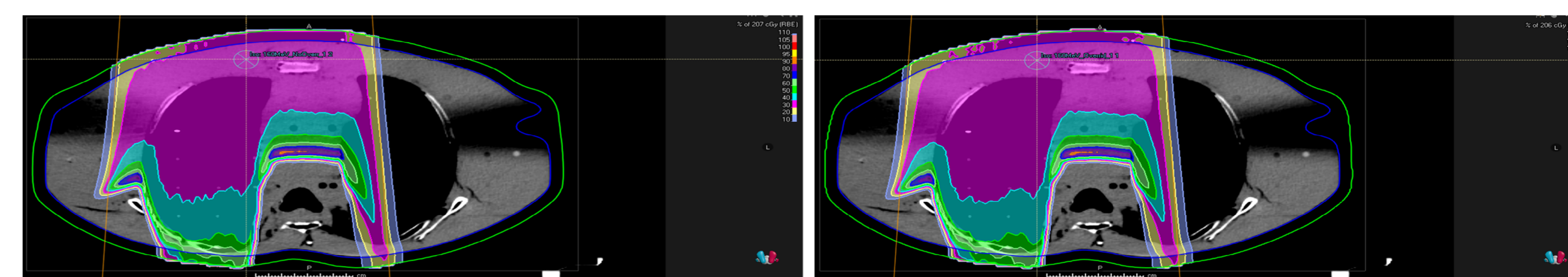


Figure 5. 160 MeV; A) no density and B) density override.

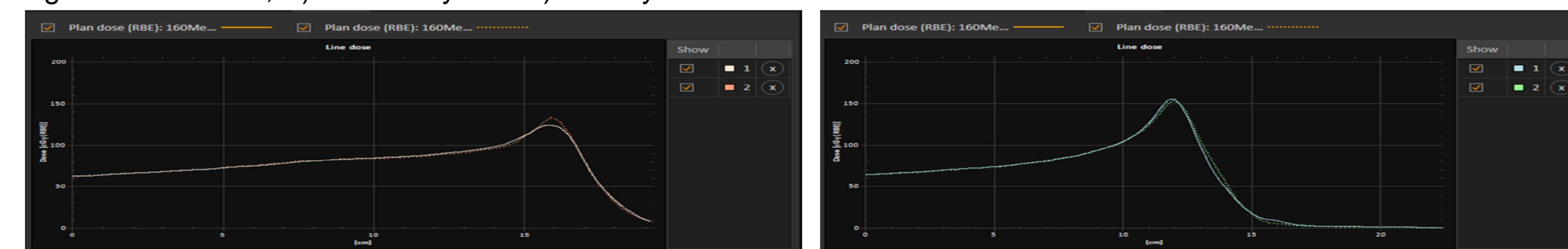


Figure 6. Dose profile perturbation at 160 MeV. A) current marker without (1) and with density override (2), B) RT-SPOT marker without (1) and with density override (2).

- To further reduce the contouring time, and keep the marker separated from the skin, prototypes with 1, 2, and 3 layers of backup material have been fabricated and compared.
- The markers were attached to an anthropomorphic phantom and scanned.
- The body contour and the markers were drawn by 3 similar expert-level dosimetrists.
- It was found that the time of contouring patient's body was reduced by 32 %.
- The time required to contour the marker wire was reduced by 23%, 47%, and 85% for markers with 1 layer, 2 layers, and 3 layers of backup material.



DISCUSSIONS:

- The average CT number for the current marker was 1378 HU (-73 to 3851 HU) and that of the new marker was 1130 HU (-27 to 2499).
- The dose perturbations due to the current marker and new marker were compared and found that the perturbations were smaller due to the new marker when the densities were not overridden to air.
- The dose perturbations were slightly larger for lower energies than those with higher energy.
- For the plan with 240 MeV energy, the dose perturbations were negligible for both the markers.
- In addition, it saves significant time in contouring the new marker.
- Also, if a planner forgets to contour and override the RT-SPOT marker, then the dose differences are smaller in comparison to the current marker.

CONCLUSION:

- These new markers are safe and reliable to use in clinics for proton therapy. This has certain advantages over the existing skin markers;
- It uses a 2.0 mm diameter, elevated, lower density, non-metallic line that minimizes artifacts.
- Existing skin markers affect dose calculations and must be manually contoured out, increasing the planning time per patient. Using this new marker saves time by reducing the need for manual contouring.
- This new marker is a flexible line that contours more easily into desired shapes and stays in place without lifting.